lumie[®]

bodyclock[®]

luxe 700FM (US)

Safety

WARNING - To reduce the risk of fire, electric shock or injury:

Read instructions for use.

If the unit is damaged, contact Lumie before use.

Keep away from water and damp. 

For indoor use only. 

The unit gets warm when the light is on, so make sure that it is not covered.

Use only the mains power adaptor supplied.

Keep the mains power adaptor and cable away from heated surfaces. If either are damaged, they should be replaced with an approved Lumie product.

This appliance does not have a on/off switch. To power down or reset, remove the mains power adaptor from the wall socket. You should therefore ensure the socket is easily accessible.

Only to be used by adults and by children when supervised.

Do not attempt to service or modify any part of the unit.

If the unit has been stored outside the recommended operating temperature range (5 °C to 35 °C / 41 °F to 95 °F), allow it to stand for at least 1 hour before applying power to avoid condensation or overheating.

Not to be used by or on individuals who are unconscious, anaesthetised or incapacitated e.g. under the influence of alcohol.

CAUTION

Do not drop the unit on the floor, hit it hard or expose it to other shocks.

Place the unit on a stable, level and non-slippery surface.

Always unplug the unit before cleaning and allow it to cool down.

Not to be used with light dimmers, timers, motion detectors, voltage transformers or extension cords.

Consult your doctor before you start using Bodyclock if you have suffered from or are suffering from serious depression.

Do not use Bodyclock as a means to reduce your hours of sleep. It is designed to help you wake up more easily, it does not diminish your need for sleep.

Intended for household use, including similar use in hotels. Do not use attachments not recommended by Lumie.

Disclaimer

It is not possible for Lumie to warn the consumer about every possible danger related to installation and use of the product. The consumer must use his or her own common sense and good judgment when installing and using the product.

	Serial number		Read instructions for use		Double-insulated (Class II electrical)
	Keep dry		For indoor use only		DC (direct current) power supply
	Manufacturer: Outside In (Cambridge) Ltd, 3 The Links, Bar Hill, Cambridge, CB23 8UD, UK.		Caution		Recyclable

Electromagnetic compatibility

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20cm between the radiator and your body. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

This Class B digital apparatus complies with Canadian ICES-003.

This device complies with Industry Canada license-exempt RSS standard(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference, and
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

For further detailed information on the electromagnetic environment in which the device can be operated refer to the website:

www.lumie.com/pages/e-m-c.

Disposal

Do not throw out the mains power adaptor or the unit with normal domestic waste. Take them to an official recycling point to be disposed of responsibly.

Contents

Introduction	6
Getting started	7
Setting the alarm time	9
Using the bedside light	11
Using the sunrise alarm	12
Using the fading sunset	14
Display brightness	15
Using the audio	16
Changing the settings	19
Sunrise options	21
Sunset options	23
Maintenance	25
Warranty	25
Technical description	26
Contact us	27

Introduction

Lumie Bodyclock wakes you naturally with a sunrise. Waking with light has been shown to improve mood, energy, productivity and the quality of sleep and awakening.¹

Wake-up lights can also help sufferers of SAD (Seasonal Affective Disorder) and winter blues cope with dark winter mornings.

At bedtime, a fading sunset helps you naturally unwind by promoting the production of the sleep hormone, melatonin.² Use together with the low-blue light feature to minimise the alerting effect of blue light and prepare for sleep.³

Lumie invented the very first wake-up light in 1993. Use all year round to encourage a regular sleep-wake pattern.

1 Acta Psychiatr Scand. 1993 Jul;88(1):67-71.


2 Sleep Biol Rhythms. 2016 Jan;14(1):47-53.

3 J Pineal Res. 2012 Aug;53(1):47-59.

Getting started

The ideal place for Lumie Bodyclock is on a bedside table near a mains socket.

- Connect the mains power adaptor to the socket in the base of the unit.
- Plug into a mains power supply and switch on.

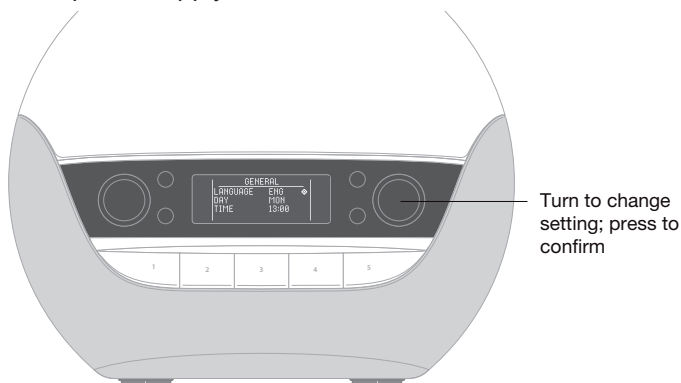The display shows the GENERAL menu.

A  cursor will prompt you to enter the essential settings:

- Turn the right dial to choose a language.
- Press to confirm and continue.

- Turn the dial to set the day.
- Press to confirm and continue.


- Turn the dial to set the hour.
- Press to confirm.
- Do the same for the minutes.
- Press to confirm and exit to the current time.

Bodyclock Luxe 700FM is now ready to use as a wake-up light, bedside light and bedside music player.

Note: The time on Bodyclock Luxe 700FM is displayed in a 24-hour format.

Front panel

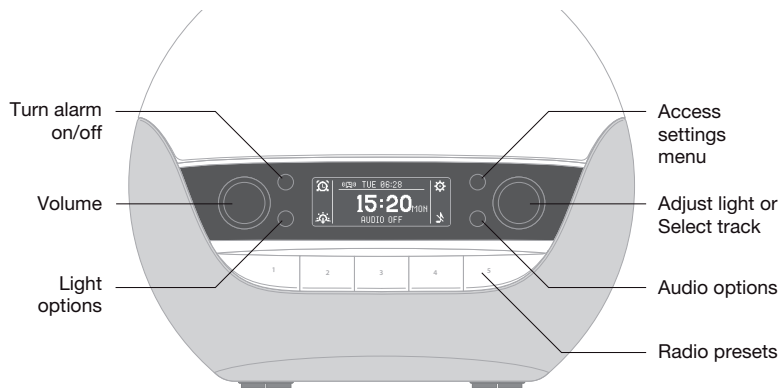The four black buttons have specific functions indicated by the icon next to them on the display. The bottom right icon changes to show the audio option playing.

The left dial adjusts the audio volume.

The right dial adjusts the light or selects audio tracks.

To adjust the light when playing audio, press the light options button  to switch to adjusting the light. To switch back to selecting tracks, press the audio options button next to the current audio icon.

The five preset buttons allow you to save and play your favourite radio stations.

Setting the alarm time

Bodyclock Luxe 700FM has two types of alarm:

With DAILY you set one alarm time and turn it on manually when you go to bed.

With WEEKLY you can set the alarm time, or choose not to have an alarm, for each day individually. For example, you may want to wake up earlier on Monday and have no alarm at the weekend.

Note: The time on Bodyclock Luxe 700FM is displayed in a 24-hour format.

To choose DAILY or WEEKLY alarm type

- Press the settings button .

The display shows the first settings menu, ALARM:

ALARM TYPE

- Press the right dial to edit ALARM TYPE.

The  cursor changes to  to show that you are editing the item.

- Turn to choose DAILY or WEEKLY and press to confirm.

The cursor will move to the next setting.

To set the DAILY alarm

If you selected DAILY use SET TIME to set an alarm time:

SET TIME

- Press the right dial to edit the hour.
- Turn to set the hour and press to confirm.
- Do the same for the minutes.
- Press the button next to  to exit from settings.

To set the WEEKLY alarm

If you selected WEEKLY you can choose for each day whether the alarm is OFF or ON, and if ON, the alarm time:

MON

- Press the right dial to edit the day; for example, MON.

The  cursor changes to  to show that you are editing the item.

- Turn to choose OFF or ON and press to confirm.

SET TIME

If you selected ON, an additional SET TIME item appears below the day.

- Turn to set the hour and press to confirm.
- Do the same for the minutes.
- Repeat for each day.
- Press the button next to  to exit from settings.

After setting the alarm time(s), the next alarm is turned on and shown on the top line of the display:

Bodyclock Luxe 700FM is preset with a 30-minute sunrise and 30-minute sunset. To use it with these settings read the following sections.

To add an audio alarm or change any other settings see *Changing the settings*, page 19, *Sunrise options*, page 21, and *Sunset options*, page 23.

To cancel an alarm

- Press the alarm button .

OFF appears at the top of the display.

If you selected WEEKLY, when the cancelled alarm time has passed the next alarm time is shown.

Using the bedside light

Lumie Bodyclock can be used as a fully dimmable bedside light.

To adjust the bedside light

Turn to adjust light;
press to switch
light colour

- Turn the right dial to adjust the brightness of the light, shown from 00 to 20 on the right of the display.

To switch to low-blue light


The blue part of the spectrum is known to have an alerting effect so it's good to minimise this at bedtime. Choose low-blue light (this looks warmer and more yellow than the white light) to optimise your bedroom lighting for relaxing sleep.

- Press the right dial to switch between normal white light and low-blue light.

Timeout: If the bedside light is not turned off manually, the light will turn off automatically after 12 hours.

Using the sunrise alarm

To wake with a sunrise alarm

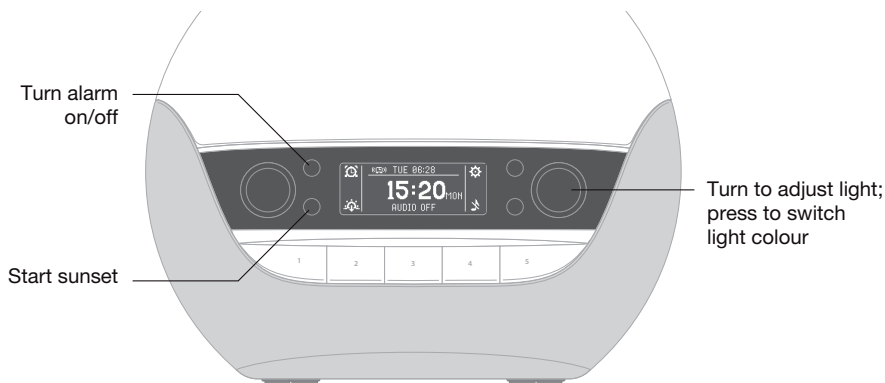If you set the alarm type to **DAILY**, you must activate the alarm before you go to sleep each time:

- Press the alarm button 

If you set the alarm type to **WEEKLY**, your alarm(s) will be turned on automatically.

The top line of the display shows the alarm time; for example:

In the morning

You'll wake with a gradual sunrise and any other sunrise options you have set. You can choose to have a wake sound or play the radio at the alarm time; see *Sunrise options*, page 21.

To snooze

- Tap the top of the lamp cover or press any preset key.

The sunrise, and audio alarm if set, will repeat after 9 minutes.

To deactivate the sunrise alarm

- Press the alarm button  or either of the dials.

This also turns off the audio alarm if you've selected one.

- Turn the right dial to turn off the light.

Timeout: If the light and audio alarm are not turned off manually, they will turn off automatically after 90 minutes.

Using the fading sunset

When you're ready to sleep, a fading sunset can help you to relax and unwind. You can choose to have the audio option you selected turn off or stay on at the end of the sunset, or gradually fade to off with the sunset; see *Sunset options*, page 23.

To start the sunset

- If you're playing audio, you may need to first press the light options button  to switch to adjusting the light.
- Turn the right dial to set the light level.
- Press the sunset button .

The display shows  to confirm the sunset is activated.

The light will start to fade and any other sunset options you have set will be activated.

To stop the sunset

- Press the sunset button .

The display shows .

Display brightness

The display automatically adjusts to match the light level in your room when the room is dark and drops to low brightness you're ready to sleep.

You can change this setting if you prefer the display to turn off completely when the room is dark:

- Press the settings button .
- Turn the left dial to select DISPLAY menu.
- Press the right dial to edit.
- Turn to set OFF or LOW.
- Press to confirm.
- Press the button next to  to exit from settings.


To check the time when the display is off

- Tap the top of the lamp cover.

The display will turn on for just a few seconds.

Using the audio

Bodyclock Luxe 700FM incorporates a high-quality audio system with FM radio that can also play music from a smartphone or other device via Bluetooth.

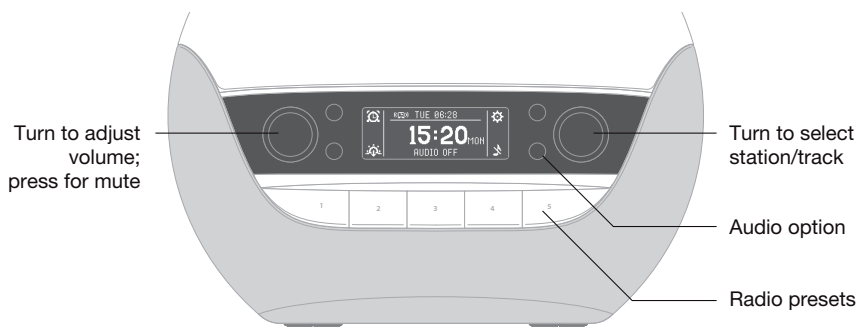To select an audio option

- If you're already listening to audio, press the audio option button to select  to turn off current audio.
- Press the audio button to choose FM radio , Bluetooth , internal sounds  or no audio .

After two seconds the audio source shown is selected.

- Use the left dial to adjust the volume, shown from 00 to 20 on the left of the display.

The following sections give more information about each of the audio options.

To play internal sounds

Bodyclock Luxe 700FM has over 20 sleep and wake sounds.

- Press the audio button until the icon next to it shows .

The display shows the SOUND menu.

- Turn the right dial to choose a sound and press to select.

- Press the button next to  to exit from the menu.

During play, you can press the right dial to go back to the SOUND menu.

To use the FM radio

- Press the audio button until the icon next to it shows **FM**.
- Turn the right dial to tune manually or press it to search for the next station.

The station or frequency is shown on the bottom line of the display:

The signal strength is shown by the  icon. If necessary adjust the telescopic aerial for best reception.

To save a station as a preset

You can save up to five stations for instant listening:

- Use the right dial to select the station you want to save.
- Press and hold a preset button until the display shows **SAVED**.
- Press the preset key to listen to your saved station.

To play the radio or music via Bluetooth

You can listen to the radio or play your own music using Bodyclock Luxe 700FM connected to your smartphone, tablet or other device via Bluetooth.

- Press the audio button until the icon next to it shows .
- On your other device, search for available Bluetooth devices and select **Lumie Luxe 700FM XXXX** for pairing (where **XXXX** refers to a random four-digit identifier).

The display will show BT PAIRED and you can now play music through Bodyclock Luxe 700FM:

Note: If you change to a different audio option, the Bluetooth pairing will be disconnected.

To turn off the audio

- Press the audio button so the icon next to it shows .


Changing the settings

This section summarises the alarm, sunrise, sunset and other options available with Bodyclock Luxe 700FM.

To access the settings menus

- Press the settings button .

The display shows the first settings menu, ALARM:

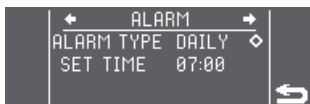- Turn the left dial to step through the different menus.
- Turn the right dial to browse the current menu.

The currently highlighted item is shown by a  cursor to its right.

Longer menus will scroll to show you all the items.

To edit a setting

- Press the right dial.

The  cursor changes to  to show that you are editing the item.

- Turn the right dial to change the item.
- Press the right dial to confirm the change.

The  cursor will move to the next item on the menu.

To exit from the settings menus

- Press the button next to .

The following table summarises each menu and menu item:

Menu	Menu items	Options	Default
ALARM	ALARM TYPE	DAILY, WEEKLY	DAILY
	MON to SUN *	OFF, ON	ON
	SET TIME	HH:MM	07:00
SUNRISE	DURATION	OFF or 15 to 90 MIN	30 MIN
	LIGHT LEVEL	1 to 20	16
	LOW BLUE	OFF, ON	OFF
	ALARM SRC (alarm source)	OFF, RADIO or SOUND	OFF
	STATION *	P1 to P5	
	TRACK *	Filename	
	VOLUME *	1 to 20	10
SUNSET	DURATION	15 to 90 MIN	30 MIN
	NIGHTLIGHT	OFF or 1 to 10	OFF
	AUDIO SRC (audio source)	OFF, MANUAL or SOUND	MANUAL
	TRACK *	Filename	
	VOLUME *	OFF, ON or FADE	OFF
DISPLAY	DISPLAY	OFF, LOW	LOW
GENERAL	LANGUAGE	ENG, FRA or ESP	ENG
	DAY	MON to SUN	MON
	TIME		13:00

* Depending on previous settings.

To change the current time


Lumie Bodyclock will keep time for 30 minutes without power. All other settings are permanently saved.

If you need to set the time use the TIME setting in the GENERAL menu.

Sunrise options


- Press the settings button .
- Turn the left dial to select SUNRISE menu.
- Continue through the following sequence to set the sunrise options:

DURATION

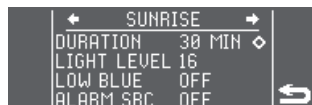- Press the right dial to edit.
- Turn to set the duration of the sunrise: OFF or between 15 and 90 minutes.
- Press to confirm and continue.


LIGHT LEVEL

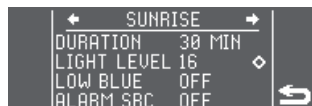Only available if you have set a DURATION.

- Press the right dial to edit.
- Turn to set the final brightness of the sunrise from 1 to 20.
- Press to confirm and continue.


LOW BLUE

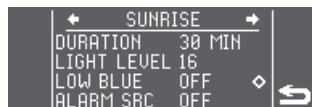Only available if you have set a DURATION.

- Press the right dial to edit.
- Turn to choose OFF for a bright white sunrise or ON for a low-blue sunrise which is warmer and more yellow.
- Press to confirm and continue.

ALARM SRC

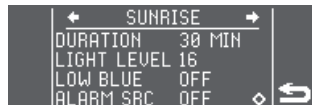- Press the right dial to edit.
- Turn to set the audio alarm option for the end of the sunrise: OFF for no audio alarm, RADIO or SOUND for one of the wake sounds.
- Press to confirm and continue.

STATION

Only available if you have set RADIO.

- Press the right dial to edit.
- Turn to set the radio preset for the end of the sunrise.
- Press to confirm and continue.

TRACK

Only available if you have set SOUND.

- Press the right dial to edit.
- Turn to set a wake sound for the end of the sunrise.
- Press to confirm and continue.

VOLUME

Only available if you have set RADIO or SOUND.

- Press the right dial to edit.
- Turn to set the volume of the audio alarm from 1 to 20.
- Press to confirm.
- Press the button next to  to exit from settings.

Sunset options

- Press the settings button .

- Turn the left dial to select SUNSET menu.

- Continue through the following sequence to set the sunset options:

DURATION

- Press the right dial to edit.
- Turn to set the duration of the sunset: between 15 and 90 minutes.
- Press to confirm and continue.

Note: The actual duration will be less than this if you start the sunset with the light already dimmed.

NIGHTLIGHT

- Press the right dial to edit.
- Turn to set a nightlight: OFF for no nightlight (sunset fades to off) or a level from 1 to 10 (sunset fades to nightlight).
- Press to confirm and continue.

AUDIO SRC

- Press the right dial to edit.
- Turn to control any audio at the start of the sunset: OFF to cancel any audio, MANUAL to leave any audio playing or SOUND to start playing one of the sleep sounds.
- Press to confirm and continue.

TRACK

Only available if you have set SOUND.

- Press the right dial to edit.
- Turn to set the sleep sound for the sunset.
- Press to confirm and continue.

VOLUME

Only available if you have set MANUAL or SOUND.

- Press the right dial to edit.
- Turn to control any audio at the end of the sunset: OFF to cancel any audio, ON to leave any audio playing or FADE to gradually decrease volume and turn off.
- Press to confirm.
- Press the button next to  to exit from settings.

Maintenance

Clean only with a dry or lightly moistened cloth. Always unplug the unit first and let it cool if necessary.

Replacement mains power adaptors are available from Lumie.

To restore the default factory settings

Only the current time will be saved.

- Turn off the power at the mains.
- Press and hold preset 3 and turn the power back on.
- Release when the display shows FACTORY DEFAULT RESTORED.

LED lifetime

The LEDs in this product will typically last for many years, even with several hours use each day. Although it's unlikely, it's possible for individual LEDs to fail. If this happens, contact Lumie; do not attempt to replace them.

Warranty

This unit is covered under warranty by Lumie LLC against any manufacturing defect for 3 years from the date of purchase. If you need to return the product for servicing, use the original packaging (including internal packing) if possible. Make sure the mains power adaptor is prevented from moving in transit and damaging the unit. If the unit is received damaged, through misuse or accidental damage, then we may have to charge to repair it; we will contact you first. Please contact Lumie for further advice.

Technical description

Subject to change without notice.

Bodyclock Luxe 700FM

Detachable accessory: Mains power adaptor (supplied) 

Input: 100–240 Vac, 50/60 Hz, 1.2 A max.; Output: 24V  1 A (24 W)

Operating conditions:

Temperature: 5 °C to 35 °C / 41 °F to 95 °F

Relative humidity: 15% to 93% (non-condensing)

Atmospheric pressure: 70 kPa to 106 kPa

Storage/transport conditions:

Temperature: -25 °C to +70 °C / -13 °F to +158 °F

Relative humidity: 15% to 93% (non-condensing)

Atmospheric pressure: 25 kPa to 106 kPa

The weights, dimensions, capacities, performance ratings, characteristics and other data displayed are mentioned only as general information. They shall not bind Lumie.

Contact us

If you have any questions or problems regarding this product, or have suggestions that you feel could help us improve our products or service please email the Lumie Careline at info@lumie.com.

Distributed in the US and Canada by:

Lumie LLC, 4600 140th Avenue N., Suite 180, Clearwater, FL 33762.

Online guide

You can view or download an online copy of the latest version of this guide on the appropriate product page of the Lumie website at <https://www.lumie.com/en-us/>.

Lumie retains the right to modify the guide as it deems appropriate.

The consumer is charged with the knowledge of the modifications set forth in the guide.

 Outside In (Cambridge) Limited, 3 The Links,
Trafalgar Way, Bar Hill, Cambridge, CB23 8UD, UK.

Lumie is a registered trademark of:

Outside In (Cambridge) Limited, incorporated in England and Wales.

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© Lumie 2023
Outside In (Cambridge) Limited, constituida en Inglaterra y Gales.